


<b>Product Name</b>	Disposable coverall		
<b>Intended use</b>	This disposable coverall is designed for Hospital, Household, Emergency, Car industry, Waste management, Gardening, Pharmaceutical, Food processing, Painting, Outing, biological chemical hazard, lab, rescue and relief, Mining, Oil and gas		
<b>Instruction for use</b>	<p>Donning</p> <p>1) Choose your correct size. Ensure that the product is suitable for the application.</p> <p>2) take out the coverall, Unfold the coverall carefully</p> <p>Check your coverall for any defects such as holes or tears and replace if damaged.</p> <p>3) Move to an area away from contamination and workplace hazards.</p> <p>Ensure you have all your personal protective equipment within easy reach.</p> <p>In a sitting position, remove footwear, bulky clothing, jewelry and anything else that may interfere with your PPE or pose a safety risk.</p> <p>4) Putting the coverall on</p> <p>Without footwear, carefully insert your leg into the coverall by pointing your foot and carefully guiding it through the ankle opening. Repeat with another leg and then pull the coverall to your waist. Carefully insert your arms to the sleeves.</p> <p>5) Roll up the legs of the coverall and put on your footwear, then roll the legs of the coverall back down to fully cover the lower leg and ankle</p> <p>Pull the hood over head, Slowly pull the zip from bottom to top and all the way up to the chin.</p> <p>Doffing</p> <p>1) Remove as much contamination as possible from your coveralls .</p> <p>2) Move to a designated area away from contamination and workplace hazards.</p> <p>Ensure that you are within easy reach of the appropriate disposal facilities.</p> <p>Ensure you have someone or something to hold on for stability.</p> <p>3) starting at the hood, pull the coverall down and step out of it.</p> <p>removing the top part of the body and arms by rolling it down and inverting the sleeves.</p> <p>roll down the coverall and step out.</p>		
<b>code</b>	OH-5000	<b>size</b>	L
<b>model</b>	Coverall with hood, double-sided tape over the entire flap. Elastic hood and ankles, cut and sewn seams, white color		
<b>Material</b>	60% Polypropylene + 40% Polyethylene, 55 ±3% gsm,		
<b>Certificate</b>	Type5 6, EN ISO 13688:2013 ,EN ISO 13982-1:2004 +A1:2010, EN 13034:2005 +A1:2009 ,EN 14325:2004 FDA, ISO12385		
<b>Packing</b>	10 pcs in a poly bag, 50 pcs/ Carton, Carton size: 50 *30*50cm Loading qty: 18000pcs in 1*20GP , 45000pcs in 1*40HQ		
<b>Shelf life</b>	3years		

		<b>Technical specification</b>					
(All measurement are for the coverall)		No.	S	M	L	XL	XX L
	chest	A	110	120	130	140	150
	Total length	B	160	165	170	175	180
	Length of sleeve	D	58	58	60	60	60
	Length of hood	G	18	18	18	18	18

orient Honest  
MEDICAL INTERNATIONAL

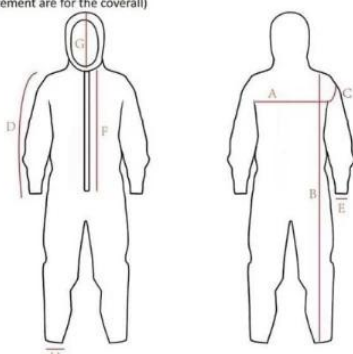

CE ISO FDA  
ISO13485

orient Honest  
MEDICAL INTERNATIONAL



CE ISO FDA  
ISO13485

 orient Honest  
MEDICAL SUPPLY CO., LTD.



**CE**  **ISO 13485** **FDA**